

PureFocus690


Hardware autofocus for Raman and IR microscopy

PureFocus690 combines advanced optics and intelligent microprocessing with a 690 nm laser. Designed primarily for infrared imaging, the system will transmit up to 1100 nm and all major Raman microscopy excitation wavelengths.

It provides a real-time focus system for infinity-corrected optical systems and offers sample detection, dynamic power adjustment, and long focus recapture range.

A motorized offset lens allows adjustment of the imaging depth into the sample, continuously holding the precise distance between the imaging focal point and a reference interface of choice or enabling different discrete image planes to be acquired using the same interface.

Designed for versatility and mounted in the infinity space, it is compatible with most brands of commercial microscopes, Prior's OpenStand imaging platform, Queensgate piezo devices, and can be adapted for third-party focus actuators.

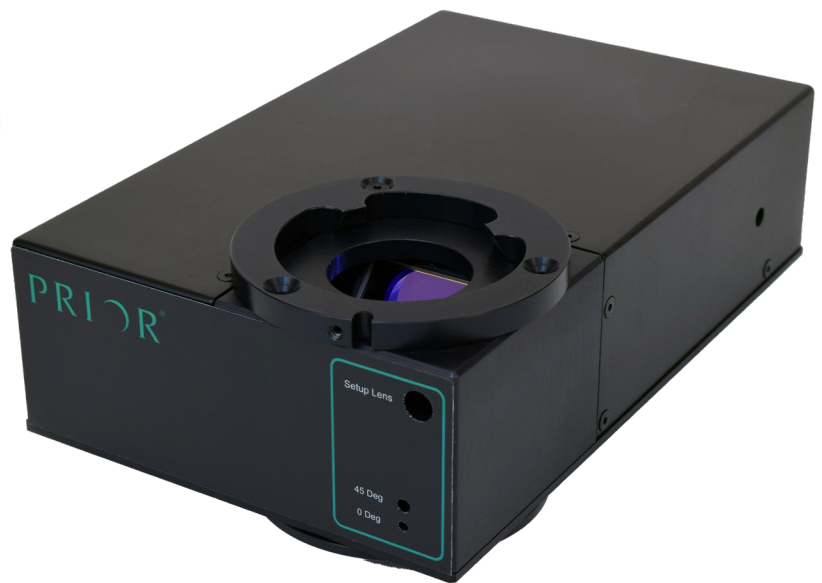



Key Features

- Rapid automated focus recovery.
- Integrated stepper motor driver and interface to Piezo and other actuators via 0-10 V analog output.
- Easily adapted to any infinity optical system, both upright and inverted.
- Stores and recalls multiple objective-specific parameters for multi-objective use.
- Customizable optics for specific imaging applications.

Applications

- Semiconductor
- Metrology
- Infrared microscopy
- Multiphoton microscopy
- Raman microscopy

Dimensions

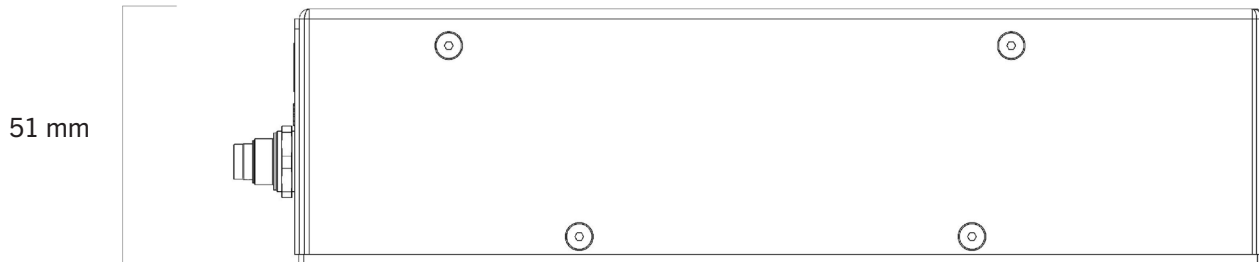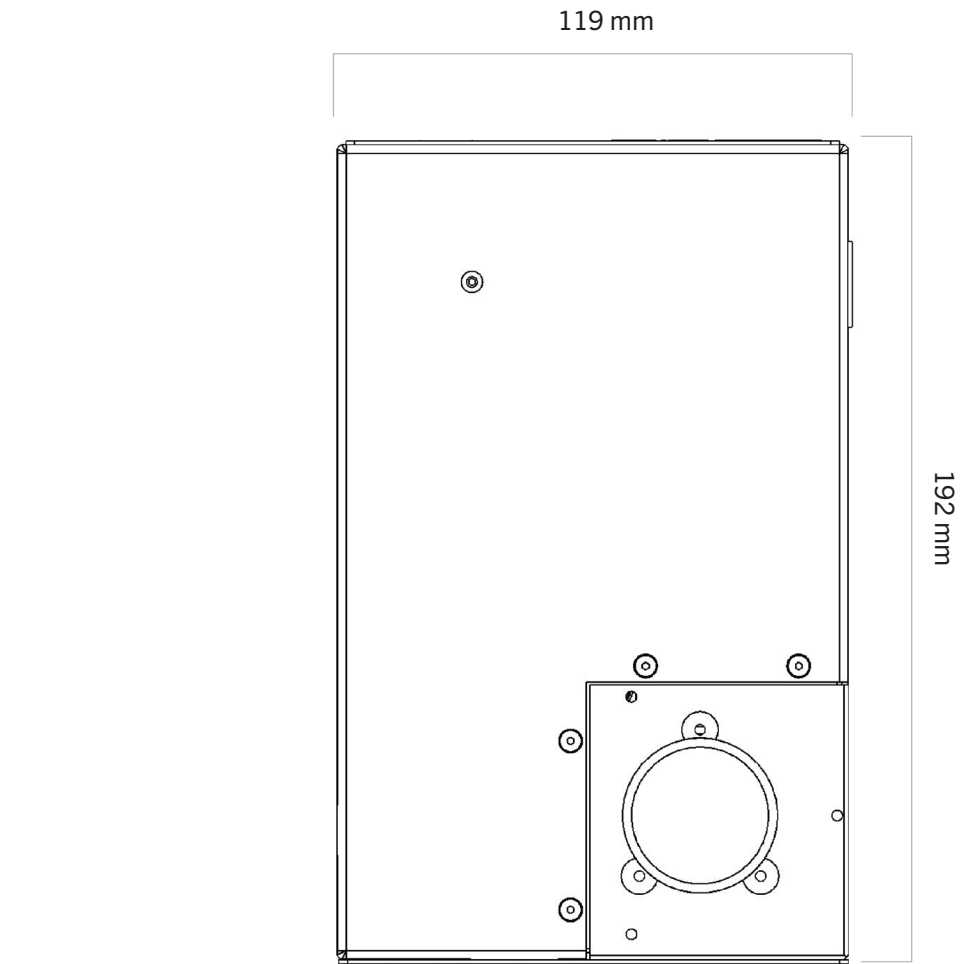

Specifications

		PF690	PF690M
Dichroic clear aperture, diameter		32 mm	32 mm
Laser optical output		< 0.77 mW	< 0.77 mW
Laser class		Class 1	Class 1
Laser wavelength		690 nm	690 nm
Focus accuracy ^{1,4}		20 % of DOF	20 % of DOF
Surface tracking method		Point on surface tracking	Averaged surface tracking
Typical update rate		1 kHz	1 kHz
Transmitted wavelength range ²		385-650 nm 730-1100 nm	385-650 nm 730-1100 nm
Capture speed ^{1,3}	5x, 0.13NA	< 30 ms	< 30 ms
	10x, 0.4NA	< 50 ms	< 50 ms
	20x, 0.5NA	< 60 ms	< 60 ms
	40x, 0.95NA	< 100 ms	< 100 ms
	100x, 0.9NA	< 1500 ms	< 1500 ms
Capture range ¹	5x, 0.13NA	+/- 15000 µm	+/- 15000 µm
	10x, 0.4NA	+/- 3000 µm	+/- 3000 µm
	20x, 0.5NA	+/- 1000 µm	+/- 1000 µm
	40x, 0.95NA	+/- 500 µm	+/- 500 µm
	100x, 0.9NA	+/- 100 µm	+/- 100 µm
Laser module mass		1.5 kg	1.5 kg

1. Representative value with optimal PID settings and signal to noise, sample and imaging system dependent.

2. >85% efficiency

3. 100 µm capture range

4. Depth of field (DOF) of the objective

Ordering Information

Part Number	Description
PF690	PureFocus690 laser module, HID controller with digipot and power supply, slice mode.
PF690M	PureFocus690 laser module, HID controller with digipot and power supply, line mode.
PF690OEM	PureFocus690 laser module, compact OEM controller and power supply, slice mode.
PF690MOEM	PureFocus690 laser module, compact OEM controller and power supply, line mode.

UNITED KINGDOM

Prior Scientific Instruments Ltd.
Units 3-4 Fielding Industrial Estate
Wilbraham Road, Fulbourn
Cambridge, CB21 5ET
United Kingdom
Email: inquiries@prior.com
Phone: +44 (0)1223 881711

U.S.A.

Prior Scientific, Inc.
80 Reservoir Park Drive
Rockland, MA. 02370
U.S.A.
Email: info@prior.com
Phone: +1 781 878 8442

GERMANY

Prior Scientific Instruments GmbH
Maria-Pawlowna-Str. 4
D-07743, Jena, Germany
Email: jena@prior.com
Phone: +49 (0)3641 242 010

JAPAN

Kayabacho 3rd Nagaoka Bldg 10F,
2-7-10, Nihonbashi Kayabacho, Chuo-Ku,
Tokyo103-0025, Japan
Email: info-japan@prior.com
Phone: +81 (0)3 5652 8831

CHINA

Prior Scientific Instruments (Suzhou) Ltd.
Room 118, Meilihua Hemu Park
No. 393 Suhong Middle Road, Suzhou Industrial Park
Suzhou, 215000, China
Email: info-china@prior.com
Phone: +86 (0)512 6617 5866



ISO 9001
Quality
Management
Systems
CERTIFIED

ISO 14001
Environmental
Management
CERTIFIED

ISO 45001
Occupational
Health and Safety
Management
CERTIFIED

